



University of
Nottingham

UK | CHINA | MALAYSIA

Atomic Force Microscopy

Force & function at the nanoscale






In this lecture...

How do we measure nanoscale forces?

Atomic Force Microscopy (AFM)

Cantilever deflection




Detection

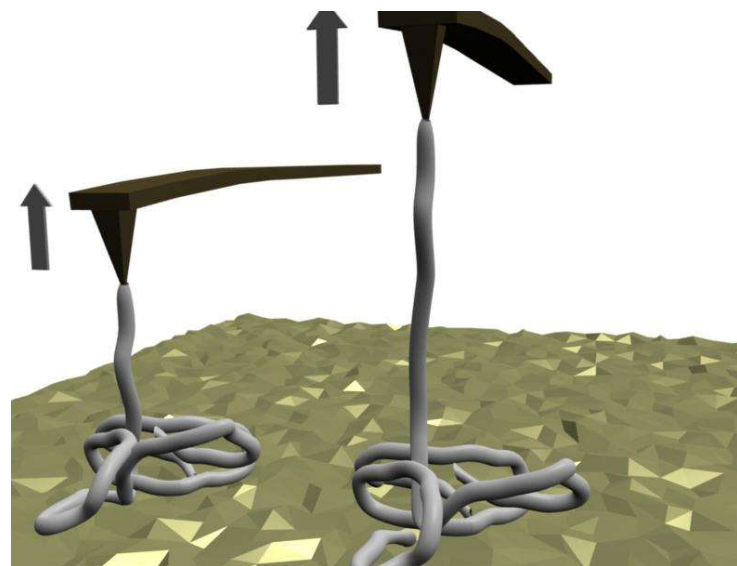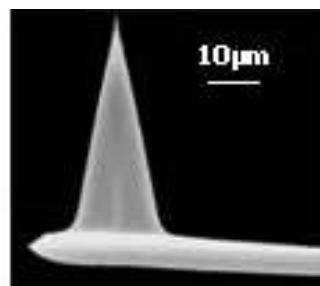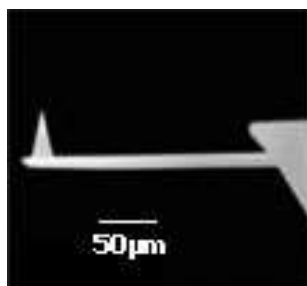



Atomic force microscopy (AFM)


Atomic force microscopy can be used to image the surface of a sample

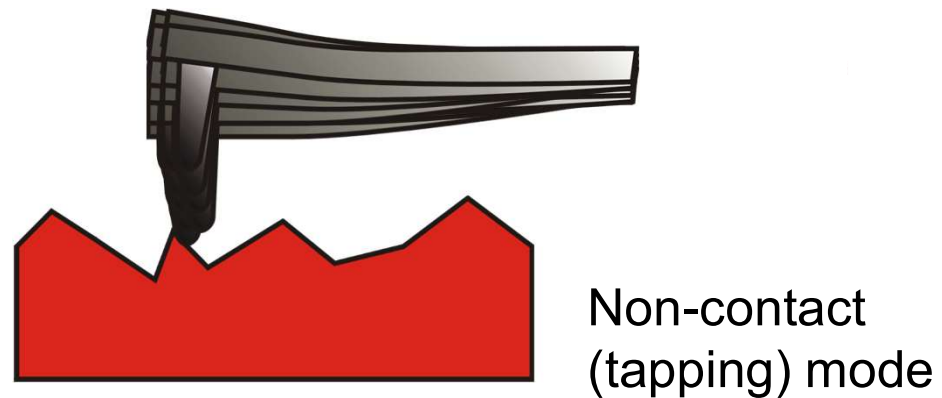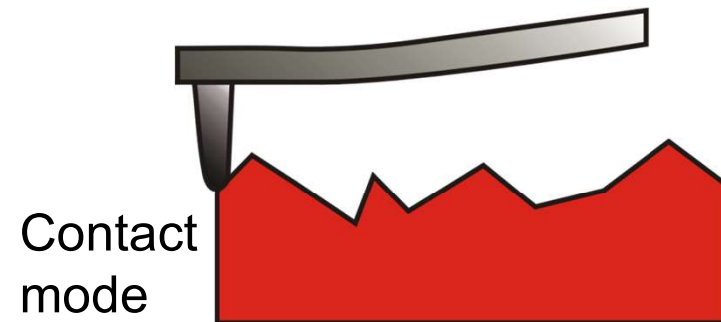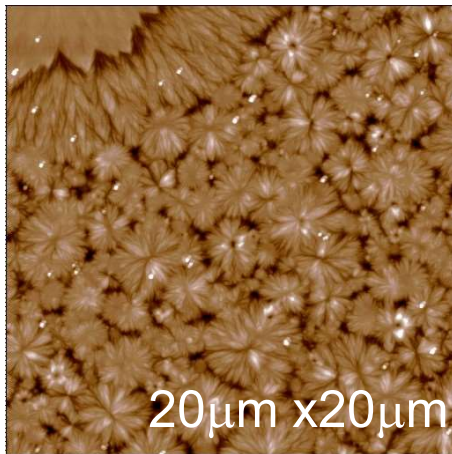



Forces at the nanoscale

Hook's law tells us we can measure a force from the stretch of a spring of known stiffness.

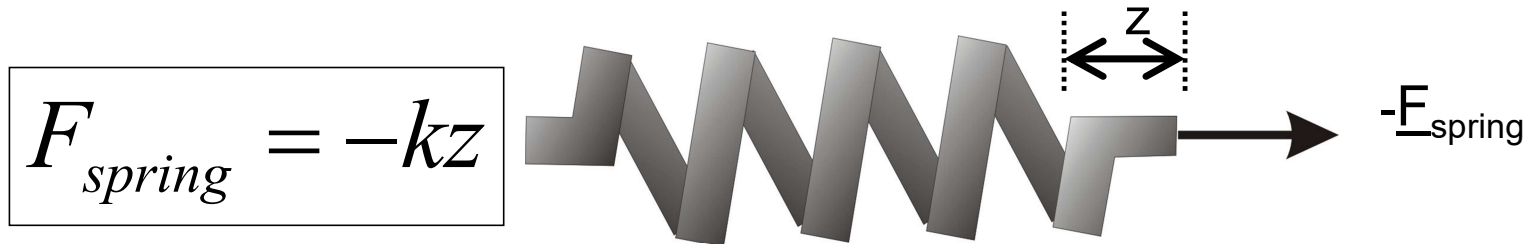


Challenges

1. How stiff is the cantilever?
2. Measure nanoscale deflections
3. Position cantilever with sub nanometre precision



Imaging bonds with Atomic Force Microscopy


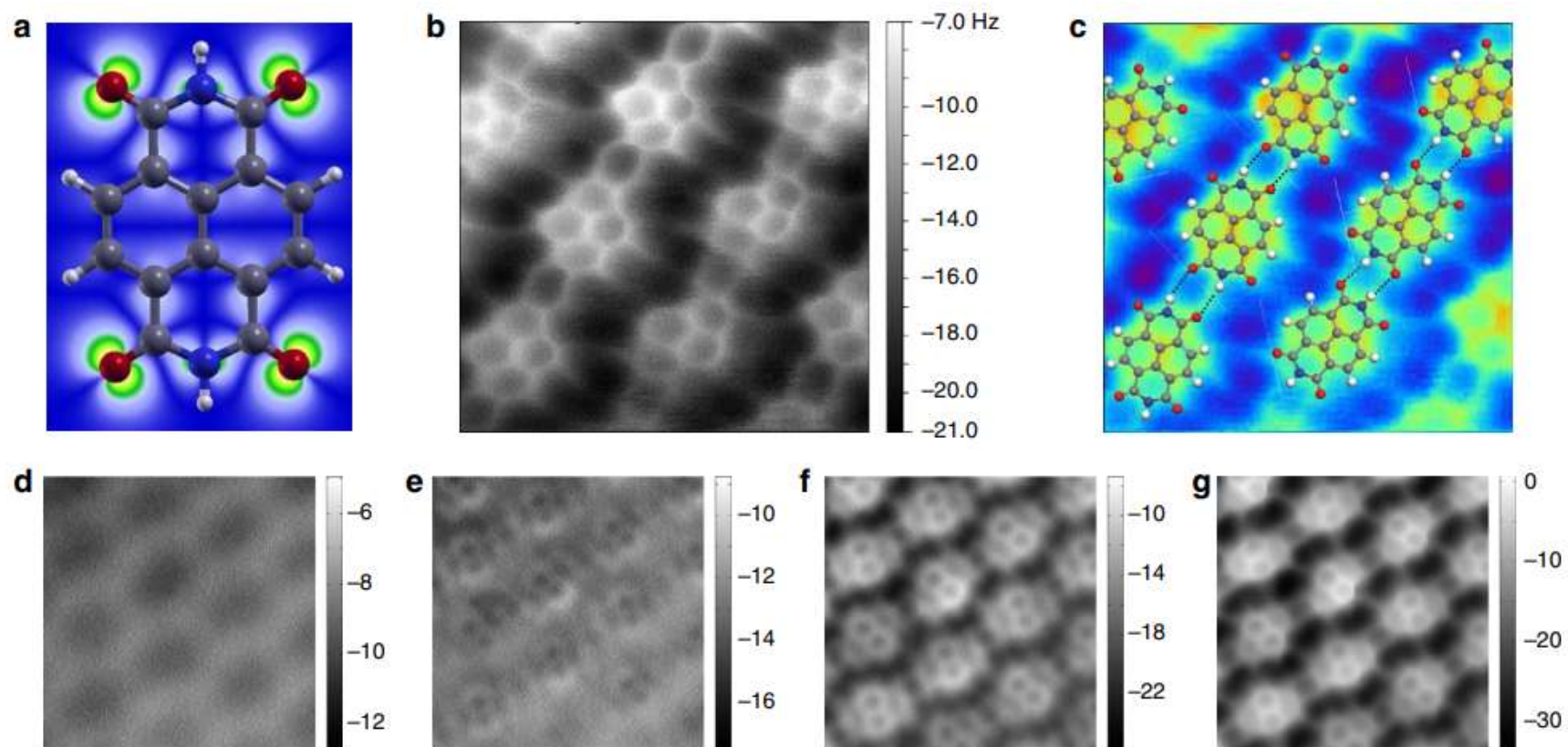


“Mapping the force field of a hydrogen-bonded assembly”
Nature Comms. 5:3931 (2014) Sweetman et al



Manipulating atoms


Switching individual silicon atoms
between two positions

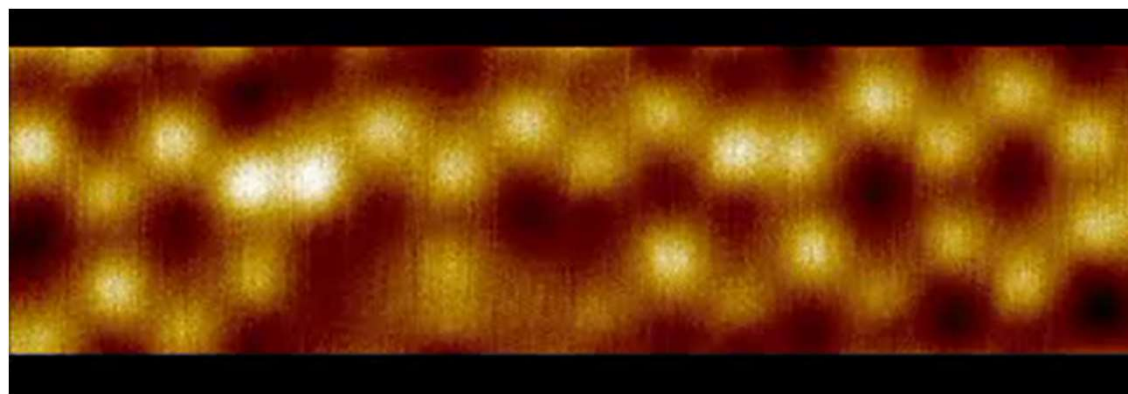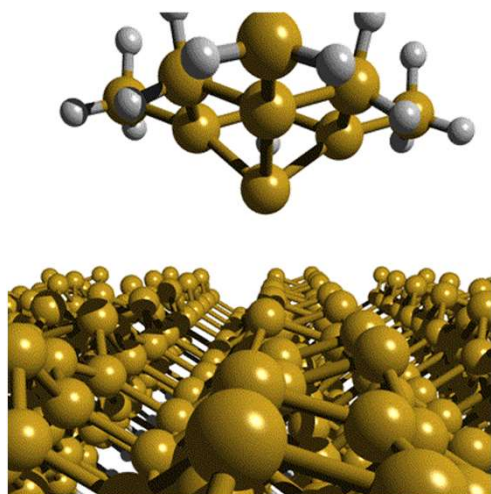



How stiff is an AFM cantilever?

The AFM cantilever bends/deflects in response to an applied force

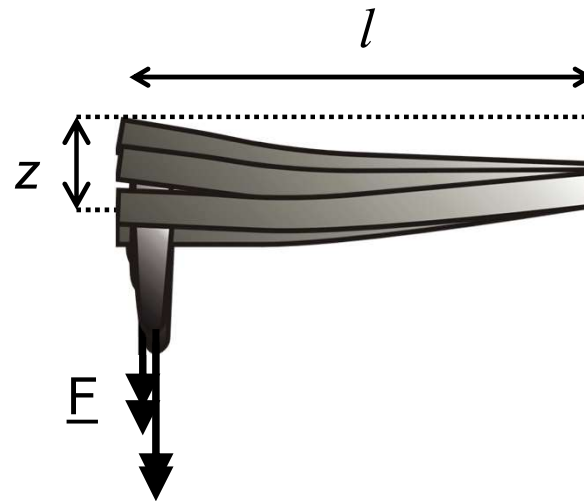


The force, F , on the cantilever is

$$F = kz$$



Derivation spring constant

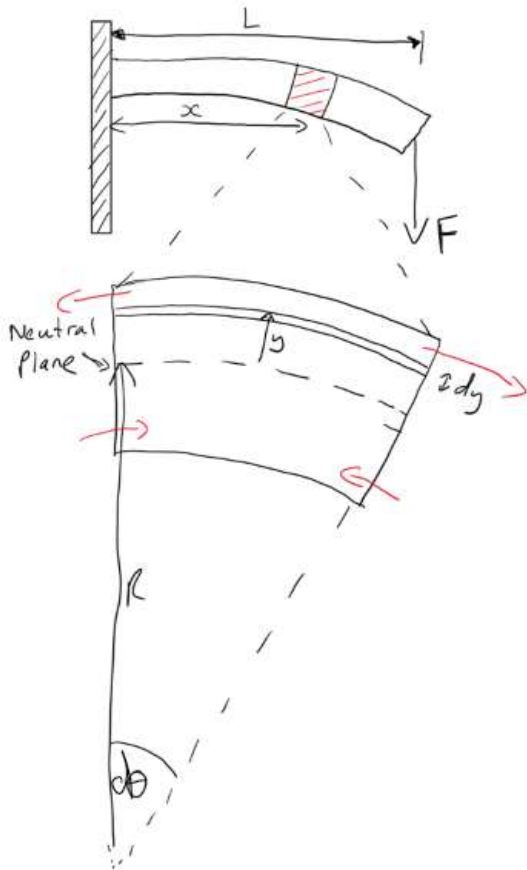

$$k = \frac{3EI}{l^3}$$

where

E = Young's Modulus of cantilever (Pa)


I = Geometric moment of inertia (m^4)



How do we detect the deflection of the cantilever?

The deflection is detected by reflecting a laser beam off the back of the cantilever.


The light is collected on a split photodiode detector

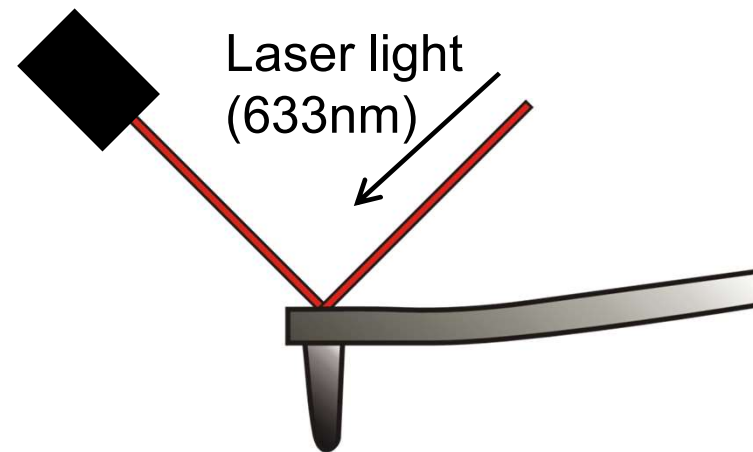



Calibrating the cantilever spring constant: thermal tuning

The spring constant of the cantilever is determined from measurements of the natural response of the cantilever to thermal fluctuations.

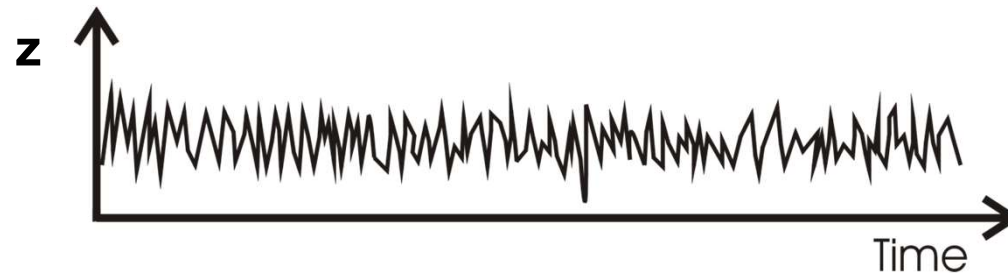

The **mean square displacement of the cantilever $\langle z^2 \rangle$** is **obtained from the noise signal** and used in combination with equipartition theory (from Thermal Physics) to give;

$$\frac{1}{2} k \langle z^2 \rangle = \frac{1}{2} k_B T \qquad k = \frac{k_B T}{\langle z^2 \rangle}$$



9.2: AFM cantilever

An AFM cantilever has a large particle stuck to its end. Far from the surface, the measured mean squared height fluctuations were 2.3nm^2 . Assuming that the tip contains most of the mass (1ng), what approximately is the resonant frequency of the cantilever?

You may assume $(2\pi f)^2 = k/m$.

$$k_B = 1.38 \times 10^{-23} \text{ JK}^{-1}, T = 300\text{K}$$





How small a force can we detect?

Typical cantilever spring constant: $0.01\text{--}50\text{ Nm}^{-1}$

We can routinely measure cantilever displacements of $\sim 0.1\text{ nm}$


We can therefore measure forces $> 10\text{ pN}$ with relative ease using AFM

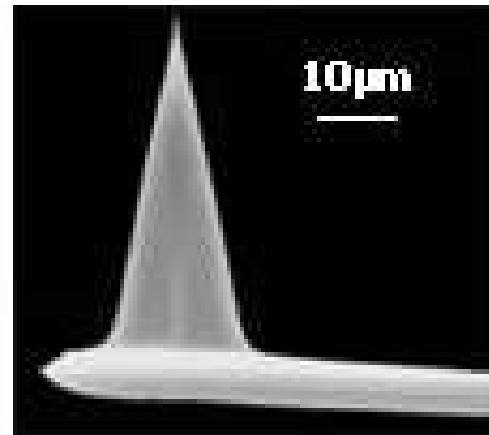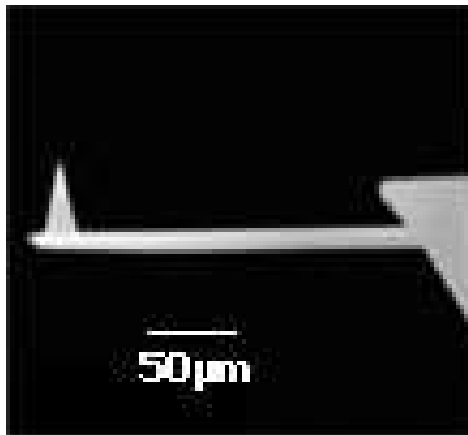



Measuring Force-Distance curves: Protein unfolding

AFM is routinely used to measure the forces associated with the folding and unfolding of individual protein molecules

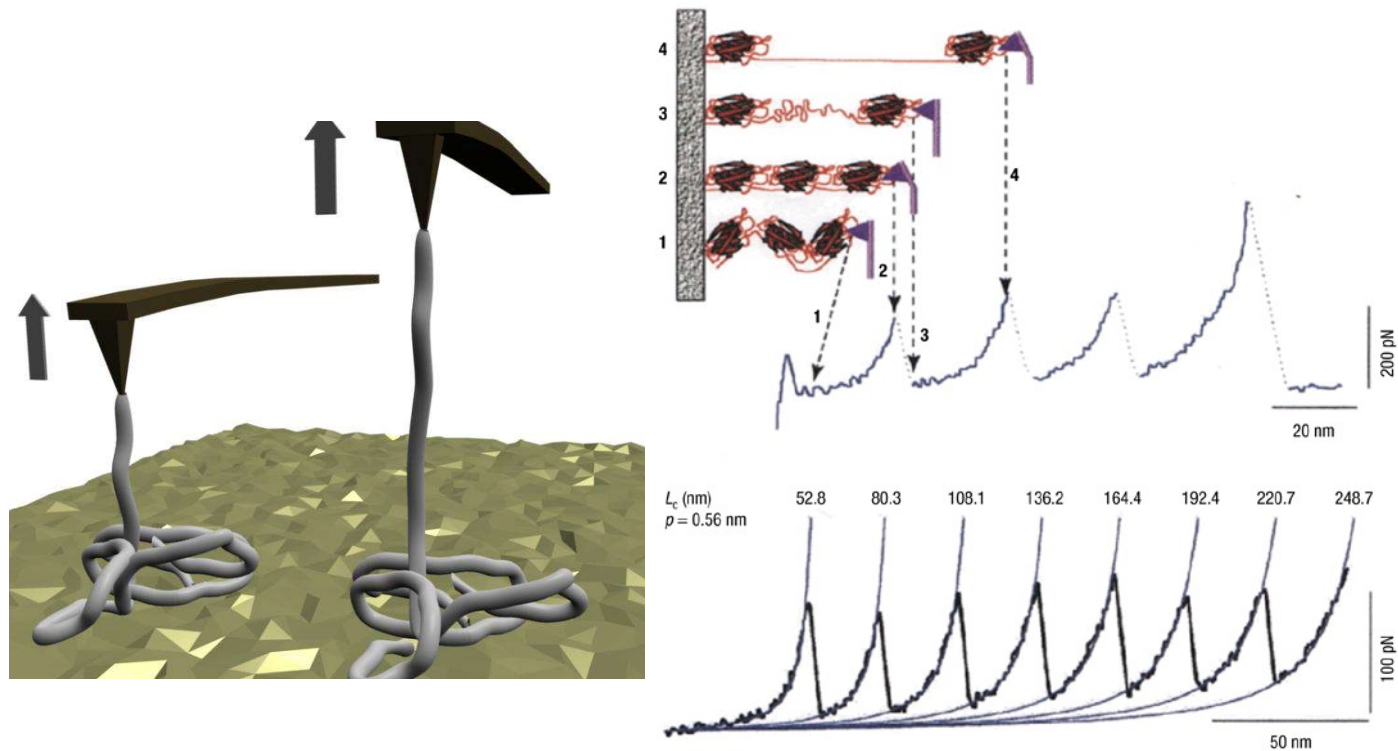


G. Bao & S. Suresh, *Nature Materials* **2**, 715 - 725 (2003)



Piezo-electric motors – nanoscale positioning of the tip

When a stress / strain is applied to a “Piezo-electric” crystal it becomes polarized producing an electric field

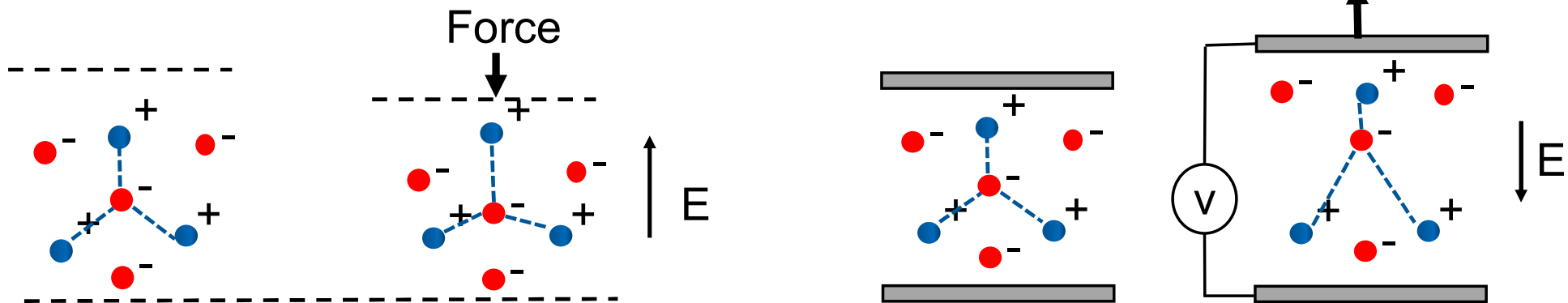

Applying an electric field can reverse this effect deforming the crystal by a small amount

By combining layers of crystal into a stack enables motion well below a nm.

$$P = -d\sigma + \alpha E$$

$$\varepsilon = \sigma/Y + dE$$




P = Polarization (dipole moment / volume), d = piezoelectric constant, σ = Stress, α = Polarizability E = Electric Field, ε = strain, Y = Young's modulus,




Mapping forces: Imaging

AFM can be used to image samples on surfaces. The cantilever is raster scanned across the surface. As the tip encounters an object the intermolecular forces deflect the cantilever.

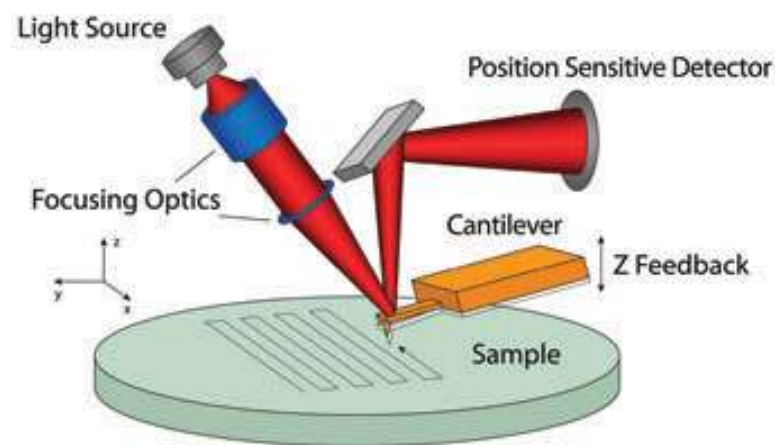

Unrestricted Optical Access from Below the Sample Plane

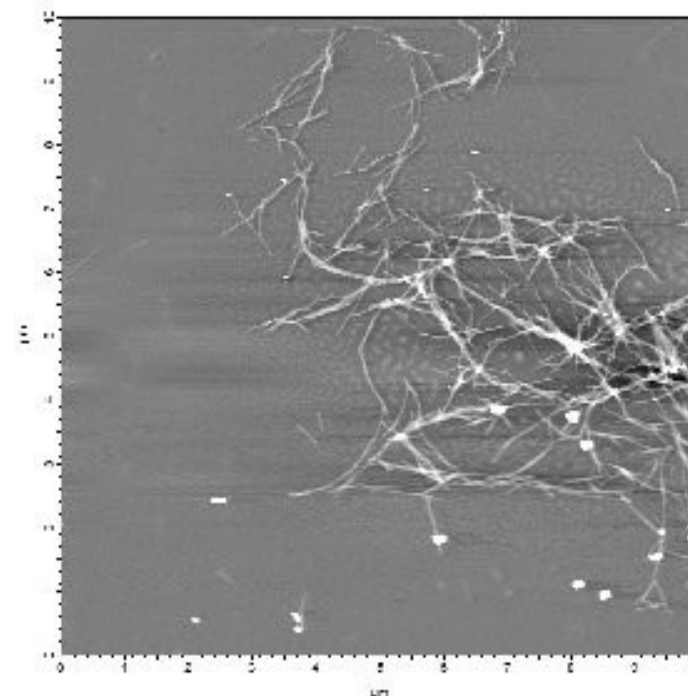


Amyloid Fibrils



Summary of key concepts

Atomic force microscopes can be used to measure forces

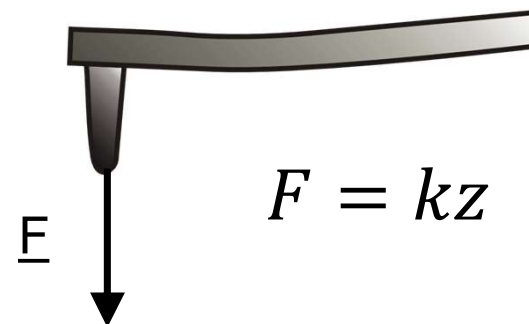

$$F = kz$$

The spring constant of an AFM cantilever is determined by its material properties and its physical dimensions

$$k = \frac{3EI}{l^3}$$

A split photodiode arrangement is used to detect deflections and to measure forces with $>10\text{pN}$ precision

